


## Jblade—Jescope video intubation laryngoscope


Jescope Ltd. recommends that all users read this manual before using the system. Failure to do so may result in injury to the patient, compromise the performance of the system, and may void the system warranty.

The Jblade slips over the connector portion of the video monitor and clicks into place, shielding the video monitor from contact with mucous membranes and non-intact skin.

Single-use Jblade are available in a variety of sizes, enabling you to treat endotracheal ranging from 5.0mmID to 8.5mmID.

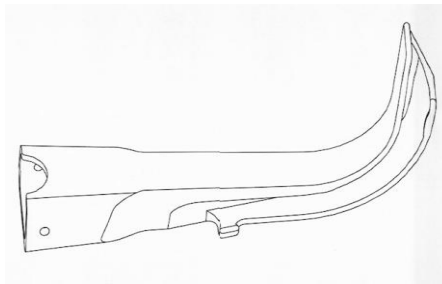



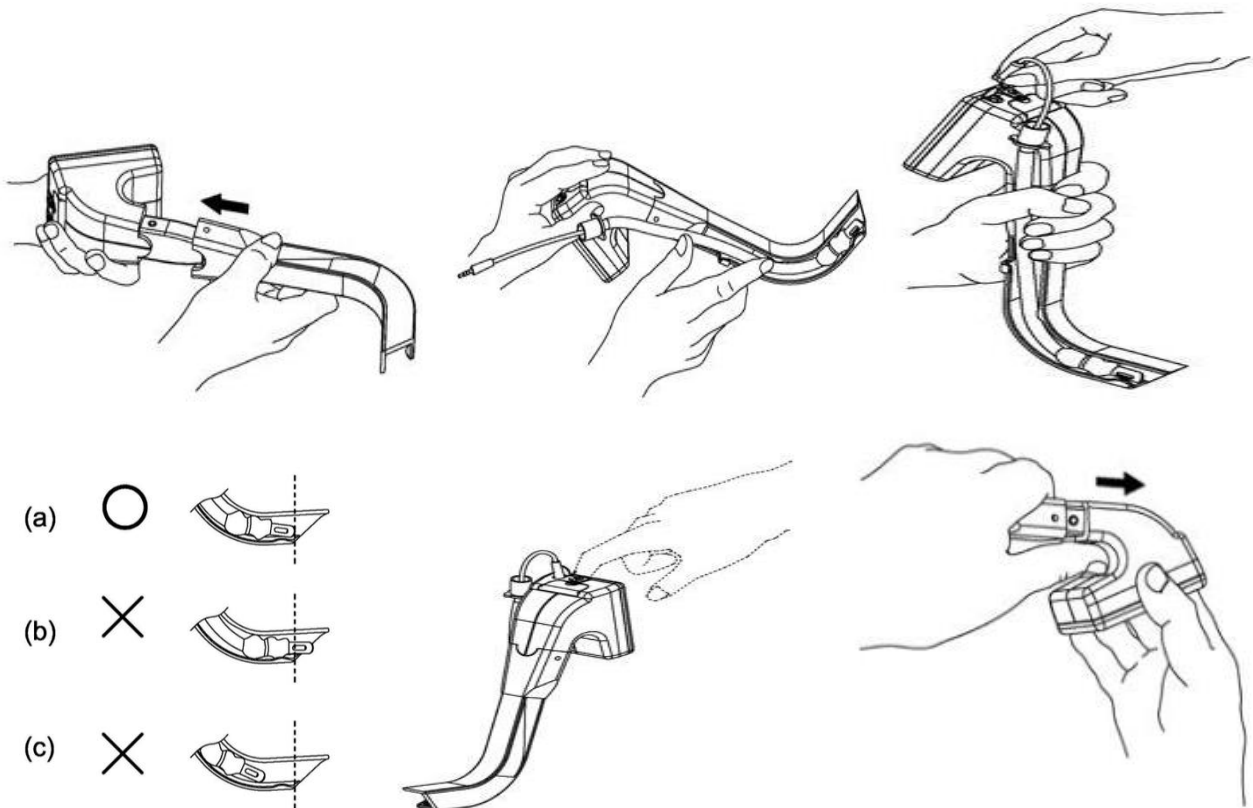









JBLADE size	Endotracheal tube size (mm I.D)
L	7.0 ~ 8.5
M	6.0 ~ 7.5
S	5.0 ~ 6.5



The Jblade is a non-sterile, single-use device. After each use, it is a biohazard, it should be removed from the video monitor and disposed of in a manner consistent with local protocols.

In addition to the user performing routine inspections before and after each use, periodic inspections should be performed to ensure safe and effective operation.

			
 J1000-BL-N	 230901		 Non-sterile
		Rx	CE
			

Report any suspected defects to Jescope Ltd. Customer Care at:

Local representative

Phone: 886-2-22984263

Email: [info@jescope.com](mailto:info@jescope.com), [www.jescope.com](http://www.jescope.com)

No.60, Dehui St., Zhongshan Dist., Taipei City 104, Taiwan (R.O.C.)